

# Highlights of the two-piece Z5s implant

## Z5s-40-12


> Implant 4 mm × 12 mm

### Screw-retained two-piece bone level implant

- Internal connection providing optimal positioning of the abutment (anti-rotational index)
- Numerous possibilities to shape the emergence profile

### Experience from over 60,000 implants placed

- Established enossal design, based on tissue-level predecessors, ensures transferable results
- Self-tapping thread with bone chip reservoir reduces insertion torque
- Excellent osseointegration through Z-Systems's proven, patented proprietary SLM® laser surface technology and gentle plasma sterilization

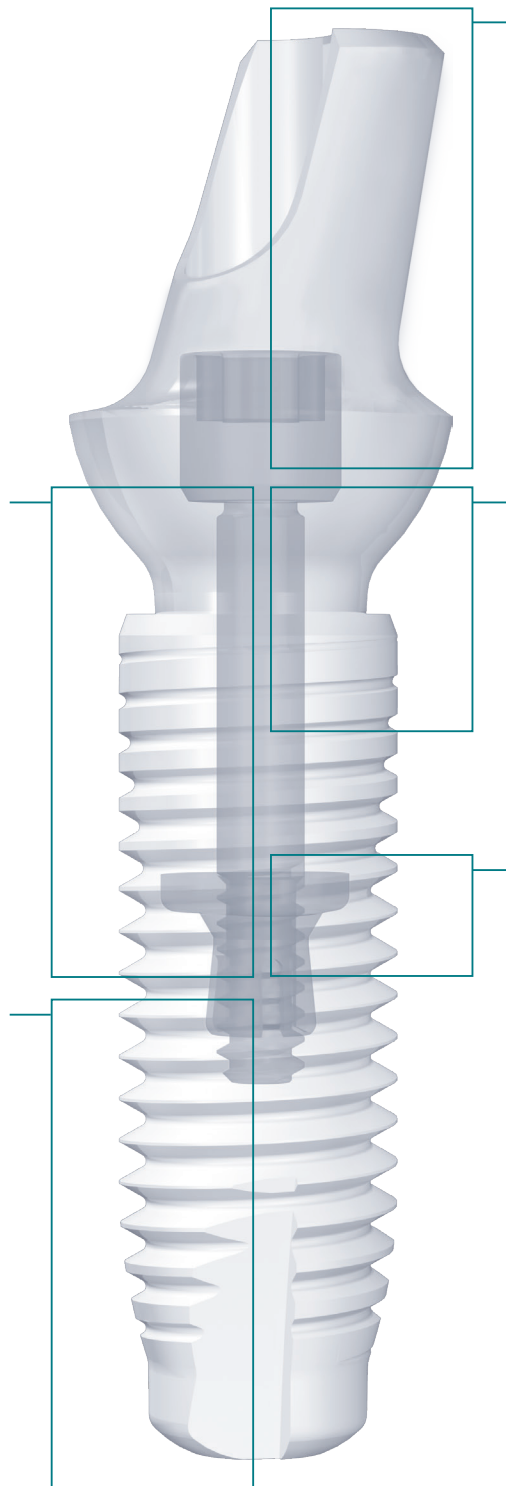

### Esthetics for professional restorations

- Straight and 15° angled abutments with gingiva height 1 mm. Other shapes coming soon - we will keep you informed.
- Variations with pre-formed mucosa profile available

### High strength connection

- Fracture resistance depending on abutment/implant diameter up to 400 Ncm (dynamic fatigue in saline solution according to DIN ISO 14801)

- Insert for secure screw retention



**Z-SYSTEMS**  
ceramic implants